

μ-ROC OEM Femtosecond Autocorrelator

ROC stands for Row Optical Correlator. Based on an ultra compact and robust inline setup, the μ -ROC takes the ROC concepts to its limit for the measurement of single-shot autocorrelation traces in the smallest housing footprint ever. Based on the most advanced innovation from Femto Easy, leveraging several years of experience in the single-shot ultrafast instrumentation, the μ -ROC is specifically designed for OEM direct integration into laser heads or laser systems.





Key features

- Ultra compact and easy to align
- Robust design, no moving parts. Non sensitive to vibrations
- Directly powered by the USB cable, no power supply required
- Suitable for any repetition rate. Single-pulse extraction possible up to 100 kHz laser repetition rate (with Trigger option)
- User-friendly and powerful software. REST API for standard software integration using simple HTTP requests

Options

Trigger


Enhanced detection

Specifications

μ-ROC Models		Ti:Sa	Ti:Sa-S	Yb	Yb-S
Pulse duration range (fs)	min	40 - 80¹	20 - 50 ¹	50 - 150 ¹	25 - 50¹
	max	1200	500	1500	1000
Accessible spectral range (nm)		700 - 900		1020 - 1080	
Input pulse repetition rate		any			
Single-pulse measurement		up to 100 kHz laser repetition rate (with Trigger option, 50 kHz without)			
Min input pulse energy ²	Single-shot	10 μJ	25 μJ	5 μJ	
	1 MHz	10 nJ	100 nJ	5 nJ	
	100 MHz	0.5 nJ	5 nJ	0.1 nJ	
Input polarization		linear horizontal or vertical			
Detection		CMOS 12 bits			
PC Interface		USB 3.1			
Beam height (mm)		20			
Dimensions (mm)		30 x 40 x 45			

² Those values give an order of magnitude. The exact sensitivity depends on many parameters (pulse duration, beam profile, wavelength...)



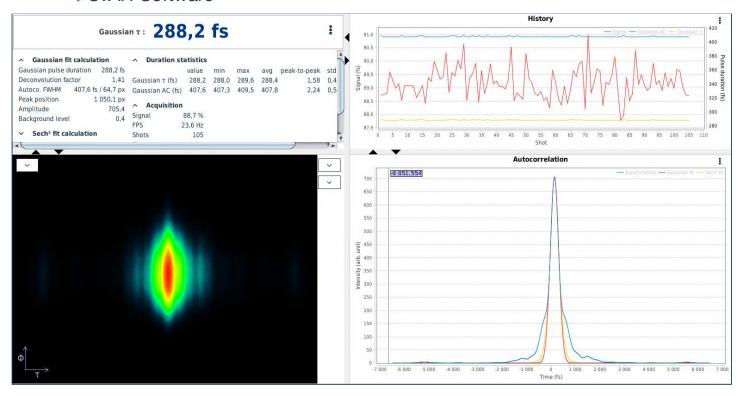

- ◆ Live extraction of shot to shot pulse duration
- ◆ Different calculation methods available for proper pulse estimation (Raw data FWHM, Gaussian fit, sech2...)
- Enhanced background & hot pixels treatment, for optimum dynamic and signal to noise ratio
- Client / Server interface, and REST API for the easiest integration
- All data exportable into most common formats



¹ The two minimum pulse duration values correspond to the Fourier limited pulse duration with and without GDD precompensation.